


Nardis ND-NX

Unison UK-Amino


150 x 4.6 mm

HPLC COLUMN APPLICATION


Aminotriazoles アミノトリアゾール類の分析

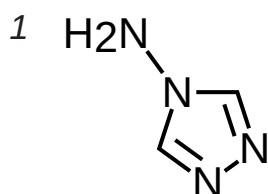


4-amino-1,2,4-triazole

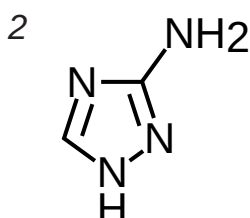

amitrole
(3-amino-1,2,4-triazole)

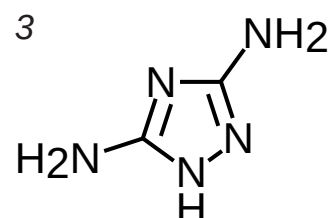

guanazole
(3,5-diamino-1,2,4-triazole)

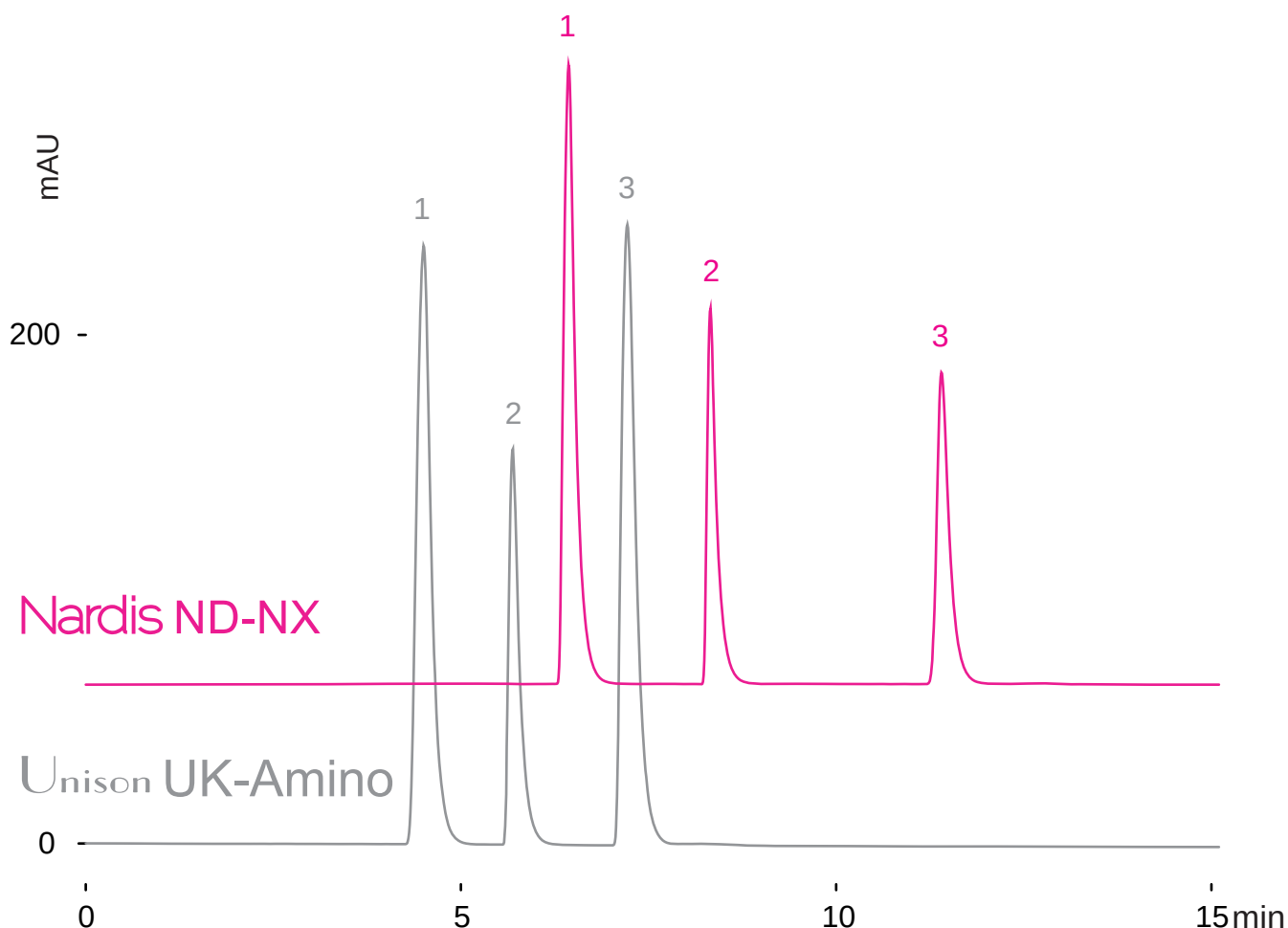

150 x 4.6 mm
A: acetonitrile
B: 100mM ammonium acetate
0 -50%B (0 -10 min), 1.0 mL/min (6 -7 MPa), 37°C
ELSD (spray chamber 70°C, drift chamber 70°C)
10 µL (3mg/mL, water)