


Unison UK-Amino

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Ethylene glycols

エチレングリコール類の LC-MS 分析

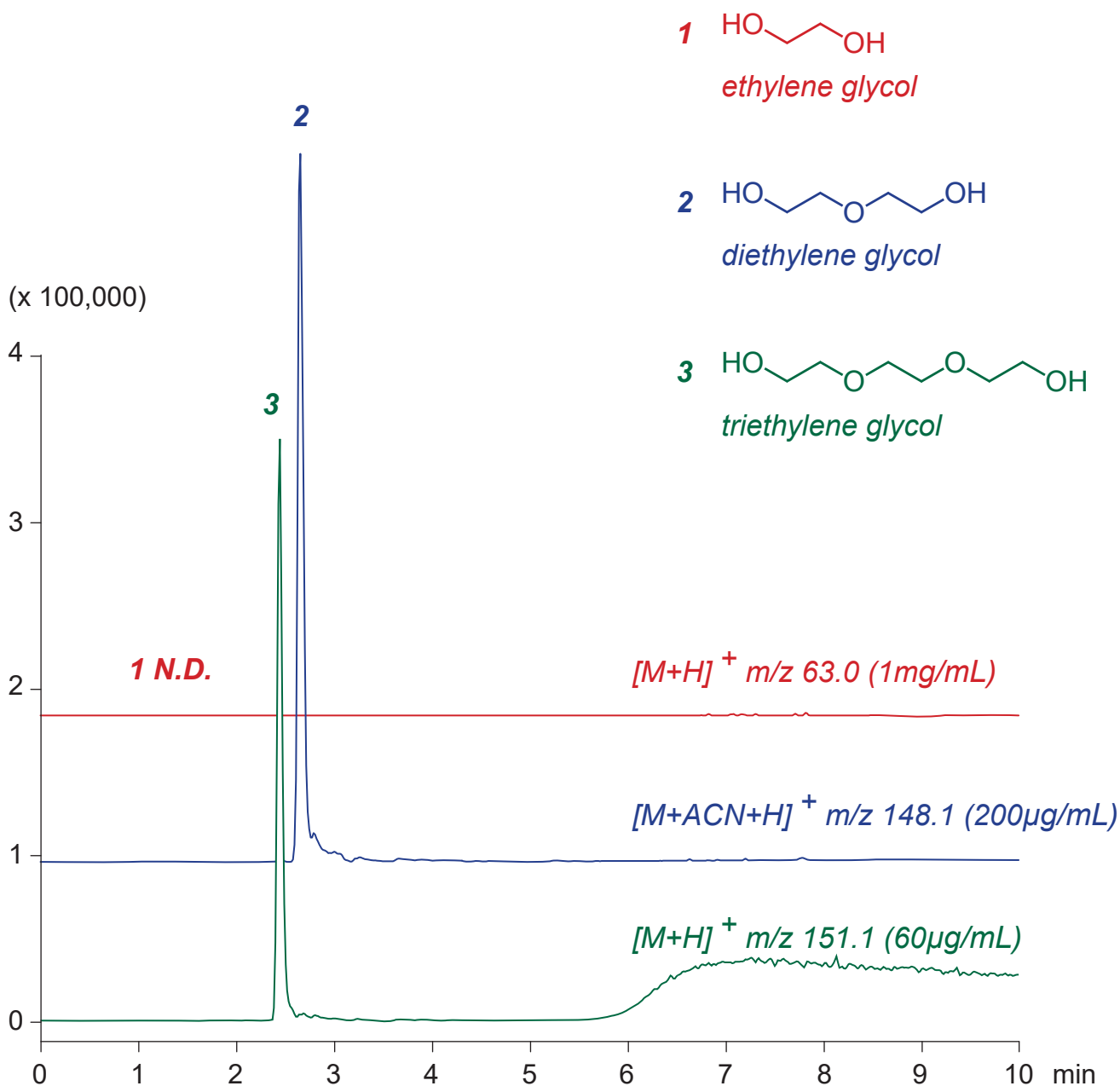

Unison UK-Amino, 150 x 3 mm

A: acetonitrile / acetic acid = 100 / 0.5

B: acetonitrile / water / acetic acid = 10 / 90 / 0.5

5%B (0-1min), 5-100%B (1-10min)

0.4mL/min (5MPa), 37°C, Single Quad. MS (ESI)

3µL (60-1000µg/mL, acetonitrile / water = 90 / 10)