


Cadenza CD-C18**Unison UK-C18**

250 x 4.6 mm

HPLC COLUMN APPLICATION**Sucrose fatty acid ester**

ショ糖脂肪酸エステル

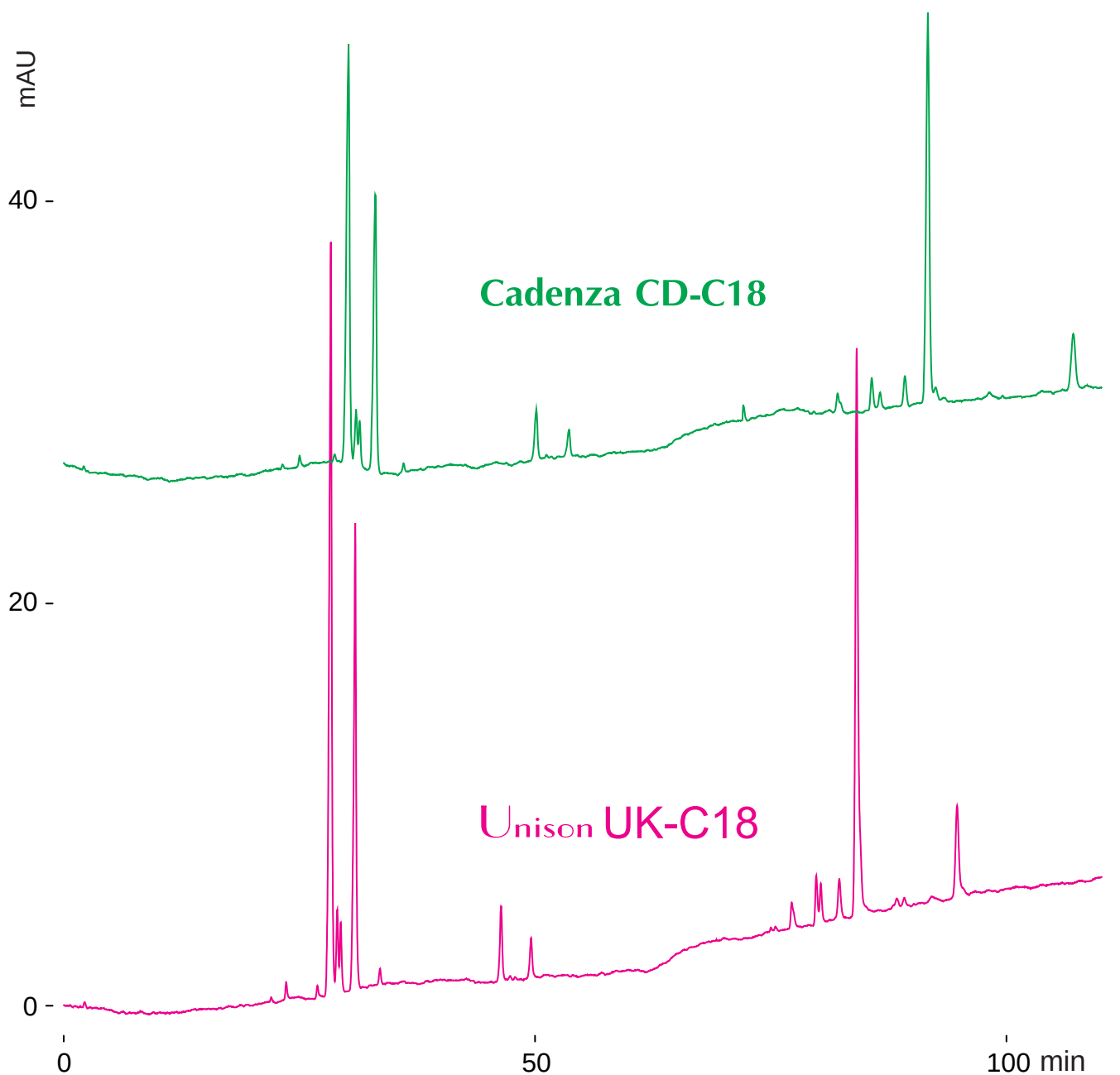

250 x 4.6 mm

A: water / acetonitrile / acetic acid = 50 / 50 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

0-40%B (0-60min), 100%B (60-110min)

1.0mL/min (17MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (10mg/mL, methanol)